

STEP 8

- Once you have selected the Therapy folder, select **Upload**.



STEP 9

- When you receive the warning box, click **Upload**.

Upload 96 files to this site:

This will upload all files from "THERAPY". Only do this if you trust the site.



STEP 10

- Click **Start** to allow data to be uploaded to your profile.

Do you confirm uploading data to Test/Test patients?



STEP 11


- Data will begin uploading. This may take a few minutes. Click **Confirm** to remove message.

The data has been uploaded and is being updated. Please check the patient data after a few minutes.



STEP 12

- After a few minutes, click **My Therapy Data** again to review your therapy data.

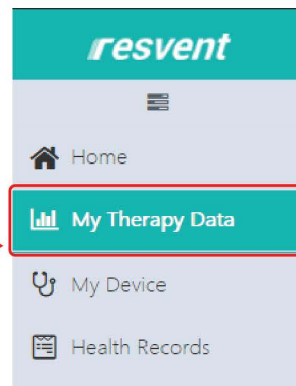

STEP 13

- You should notice that your therapy data will reflect the most current information regarding your therapy usage.
- Remove the SD Card from your PC, and insert it back into your iBreeze Device.
- You, your Physician, or Rotech staff can then create a **Custom Report**. Most insurances require an objective print out to confirm you have met compliance.
- To create a report, select the length of your report, then select **end by date**.
- Once complete, click **Create Report**.



If you have difficulty uploading your SD card, bring it to our office and we will complete the upload for you.



Failure to complete the upload will result in you being responsible to pay for the device.



*Better Sleep. Better Health.
...and it's so Easy to Do!*

Your iBreeze CPAP device records very important and helpful clinical information on a SD card - which is used by Rotech and your physician to ensure you get the best possible results and the maximum benefits from your sleep therapy.

YOUR HEALTH IS IMPORTANT!

Uploading your SD card to ResAssist each week is critical to your therapy.

Monitoring this data helps us watch for pressure inconsistencies or mask leaks & make any adjustments you may need.

Brought to you by
ROTECH[®]
HEALTHCARE INC.

LOGIN SCREEN INFORMATION

Login

Account

(Temporary password. You will change this upon login.)

Password **Sleep1**

Verification Code **348N**
↑ (EXAMPLE) ↑

☐ Remember Account

Login

Forgot Password? Register

Your machine is already registered.

Website: <https://usa.resvent.com>

My Account Username:

My Account Password:

Just follow these steps to upload your iBreeze data, and enjoy a better night's sleep




STEP 1

- Go to the ResAssist website listed above.
- Login with your Username & Temp Password.
NOTE: Your Temporary Account Password is case sensitive and will be changed when you first log in.
- Contact Rotech to obtain your account username if needed.

STEP 2

- Review the ResAssist User Services Agreement and click, **Agree**. This is done the first time you login in to the site.

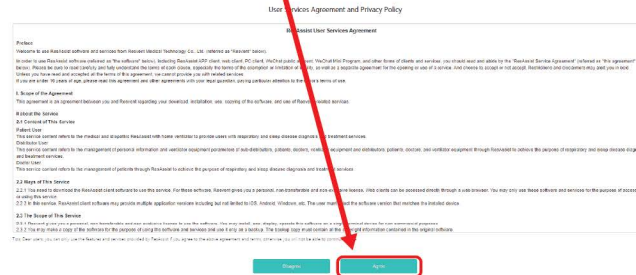

STEP 3

- On the left side of the ResAssist website, click **My Therapy Data**.



STEP 4

Remove the SD card from the iBreeze device.

- Your iBreeze system comes with an SD card already inserted in the device to store therapy information for your home care provider.

To remove the SD card:

- Stop therapy and open the device's **SD card cover**.

- Push in the SD card in to release it. Remove the SD card from the device.

Note: Do not remove the SD card during therapy.

- Do not use the SD card for any other purpose.



STEP 5

- Insert the SD card into the card reader slot on your computer.




STEP 6

- Once the SD card has been inserted into your computer, click **Upload SD Data** at the top right corner of your screen.

Upload SD Data

STEP 7

- Identify the drive for your computer's SD card reader, and click on the **Therapy** folder.

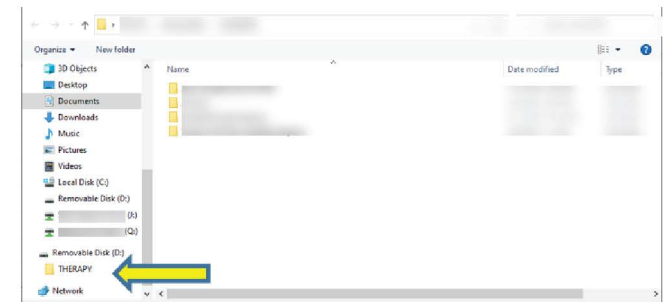

NOTE: Your computer's reader drive may be called "removable disk".