


# Step-by-Step Guide to Using the Luna QR Code Feature:


**FIRST** – Simply download and install the FREE 3B Luna QR App from the iPhone App Store or Google Play to your smart device.

*No need to upload an SD card if you follow the steps below each week!*

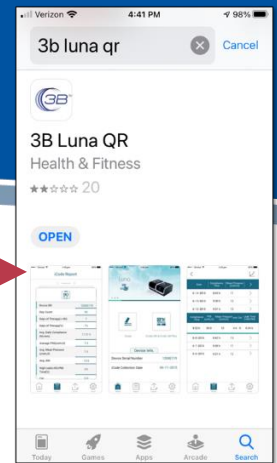

## Using the QR Code Feature – Luna II

1. Open the app and click “OK” to link the app to CPAP.
2. On the Luna II device, push alarm silence to pull up QR code. Spin dial to highlight “QR Code+”, then push the dial to select.
3. Tap iCode QR & iCode QR Plus on your phone and center QR code from device on phone.
4. System will automatically link the app to the device and upload data to the app and to iCode Connect.

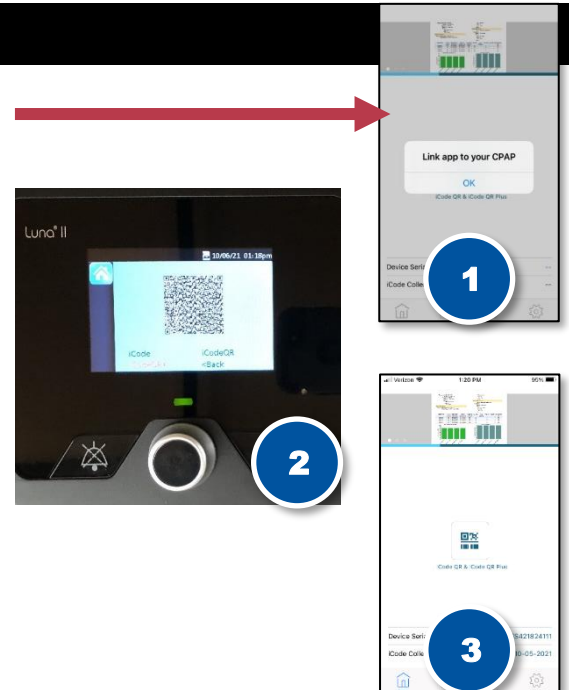

## Using the QR Code Feature - Luna G3

1. Open the app and click “OK” to link the app to CPAP.
2. On Luna G3 device, scroll to Report on the home screen, then push the dial to select.
3. Spin the dial to highlight iCode QR/QR+, then push the dial to select.
4. Tap iCode QR & iCode QR Plus on your phone and center QR code from device on phone.
5. System will automatically link the app to the device and upload data to the app and to iCode Connect.

